

Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 647

Product Details		
Size	1 mg	
Species Reactivity	Mouse	
Host/Isotype	Donkey / IgG	
Class	Polyclonal	
Туре	Secondary Antibody	
Conjugate	Alexa Fluor™ Plus 647	
Excitation/Emission Max	658/675 nm	
Immunogen	Gamma Immunoglobins Heavy and Light chains	
Form	Liquid	
Concentration	2 mg/mL	
Purification	Affinity chromatography	
Storage buffer	proprietary buffer, pH 6.5	
Contains	0.016% Methylisothiazolone, 0.016% Bromonitrodioxane	
Storage conditions	4° C, store in dark	
RRID	AB_2762830	

Applications	Tested Dilution	Publications
Western Blot (WB)	0.1-0.4 μg/mL	-
Immunocytochemistry (ICC/IF)	1-10 μg/mL	0 Publication
Miscellaneous PubMed (Misc)	-	0 Publication


Product Specific Information

To minimize cross-reactivity, the donkey anti-mouse IgG whole antibodies have been cross-adsorbed against IgG from bovine, goat, chicken, guinea pig, hamster, horse, sheep, rabbit, rat, and human. Cross-adsorption or pre-adsorption is a purification step to increase specificity of the antibody resulting in less background staining and cross-reactivity. The secondary antibody solution is passed through a column matrix containing immobilized serum proteins from potentially cross-reactive species. Only the nonspecific-binding secondary antibodies are captured in the column, and the highly specific secondaries flow through. Further passages through additional columns result in highly cross-adsorbed preparations of secondary antibody. The benefits of these extra steps are apparent in multiplexing/multicolor-staining experiments where there is potential cross-reactivity with other primary antibodies or in tissue/cell fluorescent staining experiments where there may be the presence of endogenous immunoglobulins.^M ΛM

Using conjugate solutions: Centrifuge the protein conjugate solution briefly in a microcentrifuge before use; add only the supernatant to the experiment. This step will help eliminate any protein aggregates that may have formed during storage, thereby reducing nonspecific background staining. Because staining protocols vary with application, the appropriate dilution of antibody should be determined empirically^M

Specificity: This antibody binds to heavy chains on mouse IgG and light chains on all mouse immunoglobulins. This antibody does not bind non-immunoglobulin mouse serum proteins or IgG from bovine, chicken, goat, guinea pig, hamster, horse, human, rabbit, rat or sheep.


Product Images For Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 647

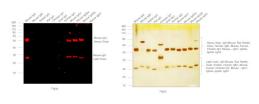

Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32787) Specificity of secondary antibody was demonstrated by specific detection of the target immunoglobulin. Antibody specificity was demonstrated by specific detection of Mouse IgG. A band at ~50 and 25 kDa corresponding to Mouse IgG Heavy and Light Chain was observed in Mouse IgG, IgG1, IgG2a, IgG2b but not in other species using Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 647 (Product # A32787) in Western Blot. {RE}


Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32787) in ICC/IF

Immunofluorescence analysis of Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 647 (Product # A32787) was performed using SH-SY5Y (positive model) and T-47D (negative model) cells stained with Nestin Monoclonal Antibody (10C2), eBioscience™ (Product # 14-9843-80). The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton™ X-100 for 10 minutes, blocked with 1% BSA for 1 hour and labeled with 2 µg/mL primary antibody for 3 hours at room temperature. Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 647 (Product # A32787, 1:2000 dilution) in 0.1% BSA in PBS for 45 minutes at room temperature, was used for detection of Nestin in the cytoskeleton (Panel a: Pink). Nuclei (Panel b: blue) were stained with Hoechst33342 (Product # H1399). F-actin was stained with Alexa Fluor® 488 Phalloidin (Product # A12379, 1:300) (Panel c: green). Panel d represents the composite image. The specificity of the secondary antibody was proved by the absence of signal in T-47D (negative model for Nestin) due to no primary antibody binding (Panel e). Nonspecific staining was not observed with secondary antibody alone (panel f). The images were captured at 40X magnification in CellInsight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) and externally deconvoluted (D.Sage et al./Methods 115 (2017) 28–41).

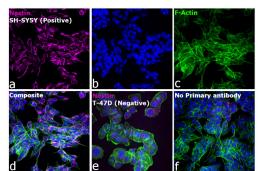

Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32787) in ICC/IF

Immunofluorescence analysis of A32790, A32816 and A32787 was performed using primary antibodies against Pericentrin (Product # PA5-53498), DDX5 (Product # PA1-31019) and ATP5A1 (Product # 43-9800) in 70% confluent log phase HEK 293 cells. The cells were fixed with 4% Paraformaldehyde, permeabilized with 0.1% Triton X-100 and blocked with 2% BSA, then incubated with the primary antibodies at 1:100 dilution each at 4 degree celsius overnight. The cells were then incubated with Donkey anti-Rabbit IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 488 (Product # A32790), Donkey anti-Goat IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 555 (Product # A32816) and Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 647 (Product # A32787) at 1:2000 dilution each in 0.1% BSA at room temperature for 45 minutes. The images were captured at 40X magnification in CellInsight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) and externally deconvoluted (D.Sage et al./Methods 115 (2017) 28-41). The specific centrosomal, nuclear and mitochondrial localization of Pericentrin (Panel a), DDX5 (Panel b) and ATP5A1 (Panel c) in the respective channels alone shows the specificity of all the 3 secondary antibodies used. Nuclei (Panel d) were stained with Hoechst33342 (Product # H1399). Panel e is the composite of Panels a-d, showing co-localisation. Panel f is control cells with

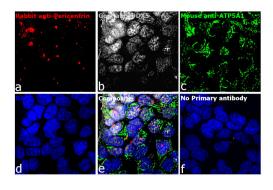

View more figures on thermofisher.cn

2

□ 170 References

Nuclear receptor-SINE B1 network modulates expanded pluripotency in blastoids and blastocysts. Nat Commun (2024)

Molecular diversity and migration of GABAergic neurons in the developing ventral midbrain. iScience (2024)

Extracellular vesicles from microglial cells activated by abnormal heparan sulfate oligosaccharides from Sanfilippo patients impair neuronal dendritic arborization. Mol Med (2024)

Cortical layer 6b mediates state-dependent changes in brain activity and effects of orexin on waking and sleep bioRxiv (2024)

Structural basis for a nucleoporin exportin complex between RanBP2, SUMO1-RanGAP1, the E2 Ubc9, Crm1 and the Ran GTPase bioRxiv (2024)

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications requised by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED. ACREMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS, SATE SLLEEPS SOLE OPTION. THERE IS NO ADMINISTRATION. THE IS NOT ADMINISTRATION. THE IS NOT ADMINISTRATION. THE IS NOT ADMINISTRATION. THE IS NOT ADMINISTRATION. TH