

Goat anti-Human IgG (H+L) Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 488

Product Details	
Size	1 mg
Species Reactivity	Human
Host/Isotype	Goat / IgG
Class	Polyclonal
Type	Secondary Antibody
Conjugate	Alexa Fluor™ 488
Excitation/Emission Max	499/520 nm
Immunogen	Gamma Immunoglobins Heavy and Light chains
Form	Liquid
Concentration	2 mg/mL
Purification	purified
Storage buffer	PBS, pH 7.5
Contains	5mM sodium azide
Storage conditions	4° C, store in dark
RRID	AB_2534080

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	0 Publication
Immunohistochemistry (Frozen) (IHC (F))	-	0 Publication
Immunocytochemistry (ICC/IF)	1-10 µg/mL	0 Publication
Flow Cytometry (Flow)	1-10 µg/mL	0 Publication
Miscellaneous PubMed (Misc)	-	0 Publication


Product Specific Information

To minimize cross-reactivity, these goat anti-human IgG (H+L) whole secondary antibodies have been affinity purified and cross-adsorbed against mouse, rabbit, and bovine serum prior to conjugation. Cross-adsorption or pre-adsorption is a purification step to increase specificity of the antibody resulting in higher sensitivity and less background staining. The secondary antibody solution is passed through a column matrix containing immobilized serum proteins from potentially cross-reactive species. Only the nonspecific-binding secondary antibodies are captured in the column, and the highly specific secondaries flow through. The benefits of this extra step are apparent in multiplexing/multicolor-staining experiments (e.g., flow cytometry) where there is potential cross-reactivity with other primary antibodies or in tissue/cell fluorescent staining experiments where there may be the presence of endogenous immunoglobulins.


Alexa Fluor dyes are among the most trusted fluorescent dyes available today. Invitrogen™ Alexa Fluor 488 dye is a bright, green-fluorescent dye with excitation ideally suited to the 488 nm laser line. For stable signal generation in imaging and flow cytometry, Alexa Fluor 488 dye is pH-insensitive over a wide molar range. Probes with high fluorescence quantum yield and high photostability allow detection of low-abundance biological structures with great sensitivity. Alexa Fluor 488 dye molecules can be attached to proteins at high molar ratios without significant self-quenching, enabling brighter conjugates and more sensitive detection. The degree of labeling for each conjugate is typically 2-8 fluorophore molecules per IgG molecule; the exact degree of labeling is indicated on the certificate of analysis for each product lot.

Using conjugate solutions: Centrifuge the protein conjugate solution briefly in a microcentrifuge before use; add only the supernatant to the experiment. This step will help eliminate any protein aggregates that may have formed during storage, thereby reducing nonspecific background staining. Because staining protocols vary with application, the appropriate dilution of antibody should be determined empirically. For the fluorophore-labeled antibodies a final concentration of 1-10 µg/mL should be satisfactory for most immunohistochemistry and flow cytometry applications.


Product will be shipped at Room Temperature.

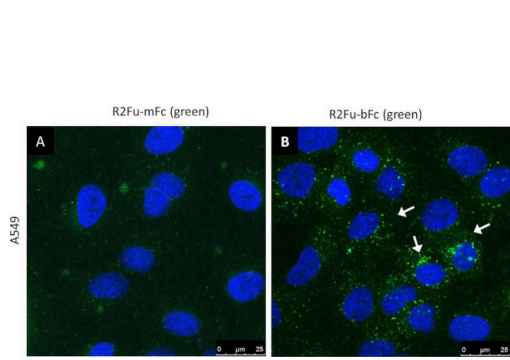

Human IgG (H+L) Cross-Adsorbed Secondary Antibody (A-11013) in ICC/IF
R2Fu-mFc and R2Fu-bFc had distinct binding profile in cells with endogenous expression of LGR4 but not in cells with endogenous expression of LGR5. (A-B), Confocal microscopy images of A549 cells incubated with R2Fu-mFc (A) or R2Fu-bFc (B). (C-D), Confocal microscopy images of SK-N-AS cells incubated with R2Fu-mFc (C) or R2Fu-bFc (D). In all images, live cells were incubated with the indicated ligands (125 nM) at 37 °C for 1 h, washed, fixed, permeabilized, and stained with Alexa488-labeled anti-human IgG. White arrows point to green punctates that represent endocytic vesicles containing RSPO ligands. Nuclei (blue) were stained with TO-PRO-3. - Image collected and cropped by CiteAb under a CC-BY license from the following publication: LGR4 and LGR5 form distinct homodimers that only LGR4 complexes with RNF43/ZNRF3 to provide high affinity binding of R-spondin ligands. <i>Sci Rep</i> (2023) Image collected and cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/37402772>), licensed under a CC BY license.

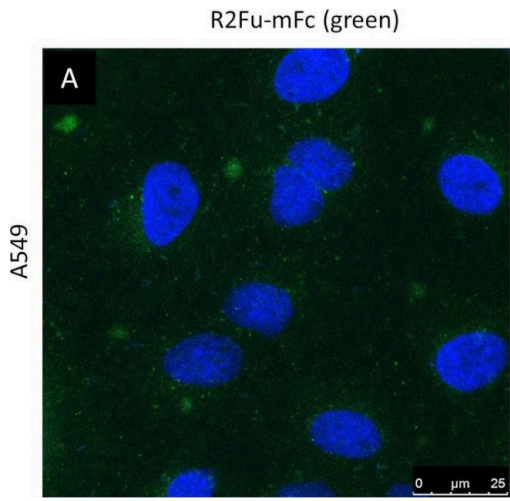

Human IgG (H+L) Cross-Adsorbed Secondary Antibody (A-11013) in ICC/IF
R2Fu-mFc and R2Fu-bFc had distinct binding profile in cells with endogenous expression of LGR4 but not in cells with endogenous expression of LGR5. (A-B), Confocal microscopy images of A549 cells incubated with R2Fu-mFc (A) or R2Fu-bFc (B). (C-D), Confocal microscopy images of SK-N-AS cells incubated with R2Fu-mFc (C) or R2Fu-bFc (D). In all images, live cells were incubated with the indicated ligands (125 nM) at 37 °C for 1 h, washed, fixed, permeabilized, and stained with Alexa488-labeled anti-human IgG. White arrows point to green punctates that represent endocytic vesicles containing RSPO ligands. Nuclei (blue) were stained with TO-PRO-3. - Image collected and cropped by CiteAb under a CC-BY license from the following publication: LGR4 and LGR5 form distinct homodimers that only LGR4 complexes with RNF43/ZNRF3 to provide high affinity binding of R-spondin ligands. <i>Sci Rep</i> (2023) Image collected and cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/37402772>), licensed under a CC BY license.

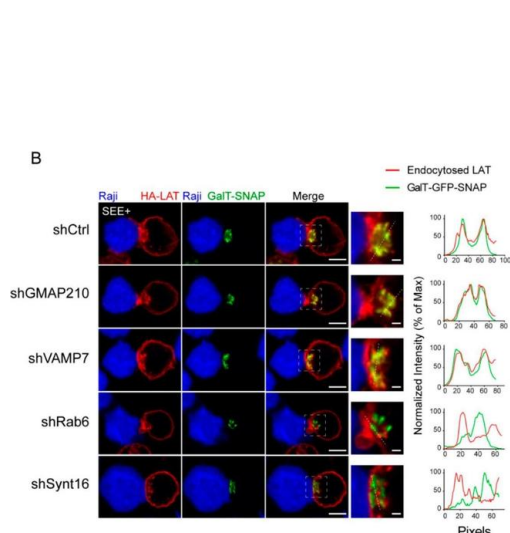

Human IgG (H+L) Cross-Adsorbed Secondary Antibody (A-11013) in ICC/IF
GMAP210 and VAMP7 do not regulate the retrograde pathway. (A) SNAP trap assay: Jurkat cells expressing GalT-GFP-SNAP and HA-LAT were stained at 4 °C with anti-HA Ab, washed, and incubated at 4 °C with BG-PEG9-NHS. After washing, cells were activated on slides for 30 min with Raji cells pulsed with SEE. Endocytosed membrane proteins, coupled with BG, covalently bound to the SNAP domain of the GalT-GFP-SNAP protein getting trapped in the Golgi. (B) Confocal images of SNAP trap assay in Jurkat cells expressing GalT-GFP-SNAP, HA-LAT and control (Ctrl) or GMAP210/VAMP7/Rab6/Synt16 specific shRNA. Labelings were performed using anti-mouse Ig (Alexa Fluor 568) to label the anti-HA Ab and anti-GFP to label the GalT-GFP-SNAP. Dashed square points to the inset localization containing the Golgi apparatus. Dashed line indicates where HA-LAT (red) and GalT-GFP-SNAP fluorescence signal intensity was measured for the plot. Images show the z-projection of summed slices from three stacks covering the Golgi apparatus. Scale bars = 5 μm. Inset scale bars = 1 μm. (C) Quantification of Mander's overlapping coefficient of GalT-GFP-SNAP over HA-LAT. (D), Mitochondrial trapping assay: Jurkat cells expressing a chimeric GFP tagged GMAP210 protein with the C-terminal hydrophobic anchor of ActA can capture vesicles interacting with GMAP210 in the mitochondria. (E) Confocal images showing the localization of VAMP7, Rab6 or Sy... Image collected and cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/33572370>), licensed under a CC BY license.

View more figures on thermofisher.cn

CD22 modulation alleviates amyloid -induced neuroinflammation. J Neuroinflammation (2025)

EphA2 regulates vascular permeability and prostate cancer metastasis via modulation of cell junction protein phosphorylation. Oncogene (2025)

A single residue switch mediates the broad neutralization of Rotaviruses. Nat Commun (2025)

Discovery of a common light chain bispecific antibody targeting PD-1 and PD-L1 by Hybridoma-to-Phage-to-Yeast (H2PtY) platform. Antib Ther (2025)

SHR-1806, a robust OX40 agonist to promote T cell-mediated antitumor immunity. Cancer Biol Ther (2024)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.