

Goat anti-Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 647

Product Details	
Size	1 mg
Species Reactivity	Guinea pig
Host/Isotype	Goat / IgG
Class	Polyclonal
Type	Secondary Antibody
Conjugate	Alexa Fluor™ 647
Excitation/Emission Max	650/671 nm
Immunogen	Gamma Immunoglobins Heavy and Light chains
Form	Liquid
Concentration	2 mg/mL
Purification	purified
Storage buffer	PBS, pH 7.5
Contains	5mM sodium azide
Storage conditions	4° C, store in dark
RRID	AB_2535867

Applications	Tested Dilution	Publications
Western Blot (WB)	1:2,500	0 Publication
Immunohistochemistry (IHC)	1-10 µg/mL	0 Publication
Immunohistochemistry (Paraffin) (IHC (P))	-	0 Publication
Immunohistochemistry (PFA fixed) (IHC (PFA))	-	0 Publication
Immunohistochemistry (Frozen) (IHC (F))	-	0 Publication
Immunohistochemistry - Free Floating (IHC (Free))	-	0 Publication
Immunocytochemistry (ICC/IF)	1-10 µg/mL	0 Publication
Miscellaneous PubMed (Misc)	-	0 Publication

Product Specific Information

To minimize cross-reactivity, these goat anti-guinea pig IgG (H+L) whole secondary antibodies have been affinity purified and cross-adsorbed against bovine, chicken, goat, hamster, human, mouse, rabbit, rat, and sheep sera prior to conjugation. Cross-adsorption or pre-adsorption is a purification step to increase specificity of the antibody resulting in higher sensitivity and less background staining. The secondary antibody solution is passed through a column matrix containing immobilized serum proteins from potentially cross-reactive species. Only the nonspecific-binding secondary antibodies are captured in the column, and the highly specific secondaries flow through. The benefits of this extra step are apparent in multiplexing/multicolor-staining experiments (e.g., flow cytometry) where there is potential cross-reactivity with other primary antibodies or in tissue/cell fluorescent staining experiments where there may be the presence of endogenous immunoglobulins.


Alexa Fluor dyes are among the most trusted fluorescent dyes available today. Invitrogen™ Alexa Fluor 647 dye is a near-infrared-fluorescent dye with excitation ideally suited to the 647 nm laser line. For stable signal generation in imaging and flow

cytometry, Alexa Fluor 647 dye is pH-insensitive over a wide molar range. Probes with high fluorescence quantum yield and high photostability allow detection of low-abundance biological structures with great sensitivity. Alexa Fluor 647 dye molecules can be attached to proteins at high molar ratios without significant self-quenching, enabling brighter conjugates and more sensitive detection. The degree of labeling for each conjugate is typically 2-8 fluorophore molecules per IgG molecule; the exact degree of labeling is indicated on the certificate of analysis for each product lot.


Using conjugate solutions: Centrifuge the protein conjugate solution briefly in a microcentrifuge before use; add only the supernatant to the experiment. This step will help eliminate any protein aggregates that may have formed during storage, thereby reducing nonspecific background staining. Because staining protocols vary with application, the appropriate dilution of antibody should be determined empirically. For the fluorophore-labeled antibodies a final concentration of 1-10 µg/mL should be satisfactory for most immunohistochemistry and flow cytometry applications.

Product will be shipped at Room Temperature.

Product Images For Goat anti-Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 647

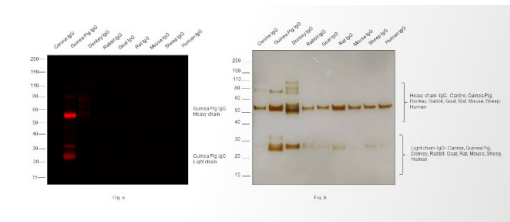


Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A-21450)
Specificity of secondary antibody was demonstrated by specific detection of the target immunoglobulin. Antibody specificity was demonstrated by specific detection of Guinea Pig IgG. Bands at ~55 and 25 kDa corresponding to Guinea Pig IgG Heavy and Light Chain were observed in Guinea Pig IgG but not in other species using Goat anti-Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 647 (Product # A-21450) in Western Blot. {RE}

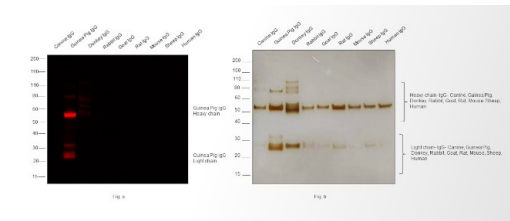

Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A-21450) in WB
Western blot was performed using Goat anti-Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 647 (Product # A-21450) and 55, 25 kDa bands corresponding to Guinea Pig IgG Heavy and Light Chain were observed in Guinea Pig but not in Canine IgG, Donkey IgG, Rabbit IgG, Goat IgG, Rat IgG, Mouse IgG, Sheep IgG, and Human IgG. Purified protein (200 ng) of Canine IgG (Lane 1), Guinea Pig IgG (Lane 2), Donkey IgG (Lane 3), Rabbit IgG (Lane 4), Goat IgG (Lane 5), Rat IgG (Lane 6), Mouse IgG (Lane 7), Sheep IgG (Lane 8) and Human IgG (Lane 9) were electrophoresed using NuPAGE™ 4-12% Bis-Tris Protein Gel (Product # NP0321BOX). Resolved proteins were then transferred onto a nitrocellulose membrane (Product # IB23001) by iBlot® 2 Dry Blotting System (Product # IB21001). The blot was probed with Goat anti-Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ 647 (Product # A-21450, 1:2500 dilution) and detected using iBright™ FL1500 (Product # A44115). Silver staining was performed to establish equivalent loading of purified proteins using the Pierce™ Silver Stain Kit (Product # 24612) (Fig b).

Guinea Pig IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A-21450) in ICC/IF
Pharmacological and genetic manipulation of KCCN channels promotes C9-ALS MNs survival and Drosophila locomotor function. A High content images of caspase 3 and MAP2 positive motor neurons from two control (MIF1 & CS14) and two C9ORF2-ALS (ALS-28 & ALS-29) lines. Scale bars: 50 μm. Bar charts represent the percentage of caspase 3 positive cells (top right) and of apoptotic cells measured by nuclear fragmentation (bottom right) in healthy and C9-ALS patient-derived motor neurons (mean ± SEM; one-way ANOVA Tukey's multiple comparison test, ****: p < 0.0001; N (replicate experiments) = 3). B High-content microscopy automated imaging quantification of caspase-3 positive cells (%) treated with DMSO or increasing concentrations of the KCNN1/3 inhibitor apamin (APA) in healthy, C9-ALS and isogenic CRISPR-Cas9 mediated deletion of the repeat expansion (ISO) patient-derived motor neurons (mean ± SEM; two-way ANOVA Tukey's multiple comparison test, ****: p < 0.0001; N (replicate experiments) = 3). C Larval crawling ability of male and female control (G4C3x3) or C9-ALS (G4C2x36) Drosophila crossed or not with two different SK mutant lines driven by nSyb-GAL4 (mean ± SEM; one-way ANOVA; ns: non-significant, *: p < 0.05, **: p < 0.01, ****: p < 0.0001; N (Drosophila larvae) > 5). D Climbing ability of male and female control (G4C3x3) or C9-ALS (G4C2x36) Drosophila crossed or not with two different SK m... Image collected and cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/34376242/>), licensed under a CC BY license.

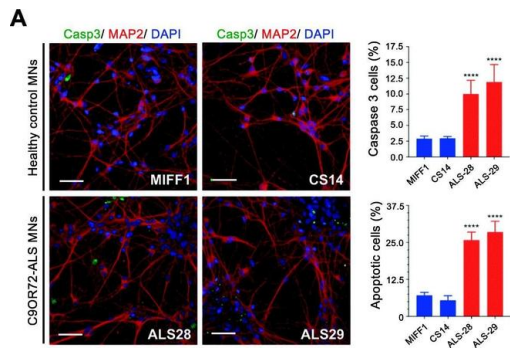

View more figures on thermofisher.cn

Intrinsic electrical activity drives small-cell lung cancer progression. Nature (2025)

Mef2c Controls Postnatal Callosal Axon Targeting by Regulating Sensitivity to Ephrin Repulsion bioRxiv (2025)

TBC1D24 interacts with the v-ATPase and regulates intraorganellar pH in neurons. iScience (2025)

Evaluation of a biomarker for amyotrophic lateral sclerosis derived from a hypomethylated DNA signature of human motor neurons. BMC Med Genomics (2025)

High resolution multi-scale profiling of embryonic germ cell-like cells derivation reveals pluripotent state transitions in humans bioRxiv (2025)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.