



# CD45RA Monoclonal Antibody (HI100), FITC, eBioscience™

<b>Product Details</b>	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG2b, kappa
Recommended Isotype Control	Mouse IgG2b kappa Isotype Control (eBMG2b), FITC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	HI100
Conjugate	FITC
Excitation/Emission Max	498/517 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11219672

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (1 μg)/test	37 Publications

#### **Product Specific Information**

Description: The HI100 monoclonal antibody reacts with human CD45RA, a 220 kDa molecule expressed by subpopulations of CD4+ peripheral T lymphocytes, CD8+ peripheral T lymphocytes, and B cells. The CD45RA+ T cell populations are mainly naive/virgin allowing the use of HI100 mAb as a phenotypic marker to discriminate T cell subsets.


Applications Reported: The HI100 antibody has been reported for use in flow cytometric analysis.


Applications Tested: This HI100 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (1 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

Excitation: 488 nm; Emission: 520 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

# Product Images For CD45RA Monoclonal Antibody (HI100), FITC, eBioscience™

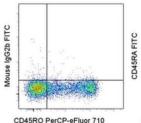

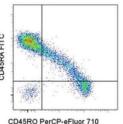

#### CD45RA Antibody (11-0458-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD45RO PerCP-eFluor® 710 (Product # 46-0457-42) and Mouse IgG2b K Isotype Control FITC (Product # 11-4732-42) (left) or Anti-Human CD45RA FITC (right). Cells in the lymphocyte gate were used for analysis.

#### View more figures on thermofisher.cn

### □ 37 References

#### Flow Cytometry (37)

Frontiers in cellular and infection microbiology

Fusobacterium nucleatum induces a tumor microenvironment with diminished adaptive immunity against colorectal cancers.

"11-0458-42 was used in Flow cytometry/Cell sorting to suggest a potential association for FN in adaptive immunity, with biological and prognostic implications."

Authors: Kim HS,Kim CG,Kim WK,Kim KA,Yoo J,Min BS,Paik S,Shin SJ,Lee H,Lee K,Kim H,Shin EC,Kim TM,Ahn JB

**Year** 2023

**Species** Human

**Nature communications** 

ADAP1 promotes latent HIV-1 reactivation by selectively tuning KRAS-ERK-AP-1 T cell signaling-transcriptional axis.

"11-0458-42 was used in Flow Cytometry to define ADAP1 as an unexpected tuner of T cell programs facilitating HIV-1 latency escape."

 $Authors: Ramirez\ NP, Lee\ J, Zheng\ Y, Li\ L, Dennis\ B, Chen\ D, Challa\ A, Planelles\ V, Westover\ KD, Alto\ NM, D'Orso\ Information (NP) and the property of the propert$ 

**Year** 2022

Species Human

View more Flow references on thermofisher.cn

## More applications with references on thermofisher.cn

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of specifications and/or accompanying package inserts ("Documentation"). No claim of studies distributed to many accompanying package inserts ("Documentation"). No claim of studies distributed to many accompanying package inserts ("Documentation"). No claim of studies distributed to many accompanying package inserts ("Documentation"). No claim of studies of the many accompanying package inserts ("Documentation"). No claim of studies distributed to many accompanying package inserts ("Documentation"). No claim of studies distributed to many accompanying package inserts ("Documentation"). No claim of studies of the warranty provided herein is valid only when used by properly trained individuals. Unless otherwise extend to any one other than the Buyer. Any mode of a sample invited to the sample invited of the sample invited invit