



Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 594

Product Details		
Size	1 mg	
Species Reactivity	Mouse	
Host/Isotype	Donkey / IgG	
Class	Polyclonal	
Туре	Secondary Antibody	
Conjugate	Alexa Fluor™ Plus 594	
Excitation/Emission Max	590/618 nm	
Immunogen	Gamma Immunoglobins Heavy and Light chains	
Form	Liquid	
Concentration	2 mg/mL	
Purification	Affinity chromatography	
Storage buffer	proprietary buffer, pH 6.5	
Contains	0.016% Methylisothiazolone, 0.016% Bromonitrodioxane	
Storage conditions	4° C, store in dark	
RRID	AB_2762826	

Applications	Tested Dilution	Publications
Western Blot (WB)	-	0 Publication
Immunohistochemistry (Frozen) (IHC (F))	-	0 Publication
Immunocytochemistry (ICC/IF)	1-10 μg/mL	0 Publication
Miscellaneous PubMed (Misc)	-	0 Publication

Product Specific Information

To minimize cross-reactivity, the donkey anti-mouse IgG whole antibodies have been cross-adsorbed against IgG from bovine, goat, chicken, guinea pig, hamster, horse, sheep, rabbit, rat, and human. Cross-adsorption or pre-adsorption is a purification step to increase specificity of the antibody resulting in less background staining and cross-reactivity. The secondary antibody solution is passed through a column matrix containing immobilized serum proteins from potentially cross-reactive species. Only the nonspecific-binding secondary antibodies are captured in the column, and the highly specific secondaries flow through. Further passages through additional columns result in highly cross-adsorbed preparations of secondary antibody. The benefits of these extra steps are apparent in multiplexing/multicolor-staining experiments where there is potential cross-reactivity with other primary antibodies or in tissue/cell fluorescent staining experiments where there may be the presence of endogenous immunoglobulins.^M


^M

Using conjugate solutions: Centrifuge the protein conjugate solution briefly in a microcentrifuge before use; add only the supernatant to the experiment. This step will help eliminate any protein aggregates that may have formed during storage, thereby reducing nonspecific background staining. Because staining protocols vary with application, the appropriate dilution of antibody should be determined empirically^M

^M

Specificity: This antibody binds to heavy chains on mouse IgG and light chains on all mouse immunoglobulins. This antibody does not bind non-immunoglobulin mouse serum proteins or IgG from bovine, chicken, goat, guinea pig, hamster, horse, human, rabbit, rat or sheep.

Product Images For Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 594

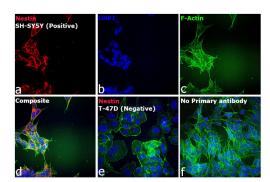


Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32744) in ICC/IF

Immunofluorescence analysis of Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 594 (Product # A32744) was performed using SH-SY5Y (positive model) and T-47D (negative model) cells stained with Nestin Monoclonal Antibody (10C2), eBioscience™ (Product # 14-9843-80). The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton™ X-100 for 10 minutes, blocked with 1% BSA for 1 hour and labeled with 2 μg/mL primary antibody for 3 hours at room temperature. Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 594 (Product # A32744, 1:2000 dilution) in 0.1% BSA in PBS for 45 minutes at room temperature, was used for detection of Nestin in the cytoskeleton (Panel a: Red). Nuclei (Panel b: blue) were stained with Hoechst33342 (Product # H1399). F-actin was stained with Alexa Fluor® 488 Phalloidin (Product # A12379, 1:300) (Panel c: green). Panel d represents the composite image. The specificity of the secondary antibody was proved by the absence of signal in T-47D (negative model for Nestin) due to no primary antibody binding (Panel e). Nonspecific staining was not observed with secondary antibody alone (panel f). The images were captured at 40X magnification in CellInsight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) and externally deconvoluted (D.Sage et al./Methods 115 (2017) 28-41).

Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32744) in ICC/IF

Immunofluorescence analysis of A32790, A32860 and A32744 was performed using primary antibodies against Pericentrin (Product # PA5-53498), DDX5 (Product # PA1-31019) and ATP5A1 (Product # 43-9800) in 70% confluent log phase HEK 293 cells. The cells were fixed with 4% Paraformaldehyde, permeabilized with 0.1% Triton X-100 and blocked with 2% BSA, then incubated with the primary antibodies at 1:100 dilution each at 4 degree celsius overnight. The cells were then incubated with Donkey anti-Rabbit IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 488 (Product # A32790), Donkey anti-Goat IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 680 (Product # A32860) and Donkey anti-Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor™ Plus 594 (Product # A32744) at 1:2000 dilution each in 0.1% BSA at room temperature for 45 minutes. The images were captured at 40X magnification in CellInsight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) and externally deconvoluted (D.Sage et al./Methods 115 (2017) 28-41). The specific centrosomal, nuclear and mitochondrial localization of Pericentrin (Panel a). DDX5 (Panel b) and ATP5A1 (Panel c) in the respective channels alone shows the specificity of all the 3 secondary antibodies used. Nuclei (Panel d) were stained with Hoechst33342 (Product # H1399). Panel e is the composite of Panels a-d, showing co-localisation. Panel f is control cells with

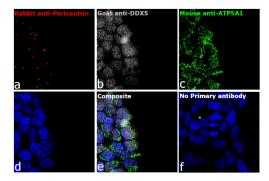

Mouse IgG (H+L) Highly Cross-Adsorbed Secondary Antibody (A32744) in ICC/IF

Immunofluorescent analysis of tubulin in U2OS cells. The cells were fixed with 4% formaldehyde for 20 mins, permeabilized with 0.5% Triton X-100 in PBS for 20 mins, washed 3X in PBS and blocked with 3% BSA in PBS for 30 mins at RT. Cells were stained with a tubulin antibody at a dilution of 1:2000 in 3% BSA in PBS for 1 hr at RT, washed 3X in PBS and then incubated with Invitrogen Alexa Fluor Plus 594 donkey anti-mouse IgG secondary antibody (Product # A32744) prepared in 3% BSA in PBS at a dilution of 1:1000 for 1 hr at RT in the presence of NucBlue Live ReadyProbes Reagent (Product # R37605). The image contains overlay of tubulin (red) and nuclei (blue). Images were taken on an EVOS FL Auto 2 Imaging System (Product # AMAFD2000) with an Olympus 40X Super Apochromat objective (Product # AMEP4754) at 40X magnification. Actin was stained using Alexa Fluor Plus Phalloidin (Product # A30105)

View more figures on thermofisher.cn

□ 145 References

Limitations of a proper SFTSV mouse model using human C-type lectin receptors. Front Microbiol (2025)

Neuronal guidance factor Sema3A inhibits neurite ingrowth and prevents chondrocyte hypertrophy in the degeneration of knee cartilage in mice, monkeys and humans. Bone Res (2025)

Suppression of ADP-ribosylation reversal triggers cell vulnerability to alkylating agents. Neoplasia (2025)

Loss of Dnah5 Downregulates Dync1h1 Expression, Causing Cortical Development Disorders and Congenital Hydrocephalus. Cells (2024)

Mice lacking acid-sensing ion channel 2 in the medial prefrontal cortex exhibit social dominance. Sci Adv (2024)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications required by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, RTINESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERFIOIS DISLIMITED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, RTINESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE EXEMBLY FOR NON-CONFORMING PRODUCTS OR NON INFRICE IS NO BELIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS; SA THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or any type of consumption by